

## Nuance Universal & Flowable Composites

Nuance Universal is a light-cured micro hybrid composite that simplifies the composite restoration process like never before. It's translucent when thin and opaque when thick. One body shade can be used for an entire restoration with flawless, natural-looking results. The multi-faceted fillers reflect light from the surrounding tooth structure and create a seamless integration, making restorations nearly invisible. With excellent physical properties like 1.9% volumetric shrinkage and state-of-the-art compressive strength, flexural strength and modulus, Nuance is a true universal that works for all classes of restorations.



- Superb handling properties
- Creamy consistency easily spreads
- Available in nine shades
- Convenient flip-top syringe or pre-loaded tips for minimal waste




Before



After

### Volumetric Shrinkage (%)\*

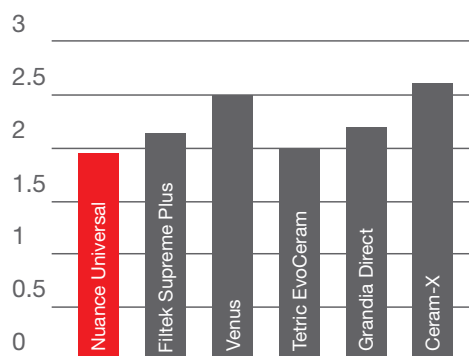

\*Data on file.

### Nuance FLOW

Nuance is also available in a highly polishable flowable composite.

- Easy to place without running or dripping
- Unique self-leveling viscosity
- Radiopaque
- Convenient syringe delivery
- Available in eight shades
- Ideal for the pulpal or axial layer in posterior restorations



Available through

 **HENRY SCHEIN®**  
DENTAL

To Order: 1-800-372-4346 8:00am–8:30pm (et)  
[www.henryscheindental.com](http://www.henryscheindental.com)

# nuance®

# nuance® FLOW



## Features

Translucent when thin and opaque when thick

Contains multi-faceted fillers

Micro hybrid composite

Syringe delivery system

## Benefits

Single shade can be used for the entire restoration with flawless, natural looking results

Fillers reflect and refract light, mimicking the opacity, translucency and optical characteristics of tooth structure

Delivers great esthetics

Innovative design makes dispensing and working with composite material easier and more intuitive

### Hinged Cap

Designed with an attached, flip-top cap that can't become lost or separated from the syringe.



### Stands on Your Tray

Nuance easily stands on your bench top for easy visibility and accessibility.



### Color-Coded


Shade Identification

Each Nuance syringe has a clearly labeled shading color-code, making it fast and easy to find on the shelf or in the caddy.




### Audible and Tactile Feedback for Unit-Dose Application

With every half-rotation of the Nuance syringe, the clinician will hear an audible click confirming that 0.1 grams of composite has been dispensed.

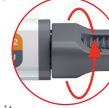


### Easy-to-See Tick Marks

The clinician or dental assistant can tell at a glance how much Nuance composite is remaining in the syringe so no one is caught short during a procedure. Experience just-in-time reordering to keep exact inventory on hand.

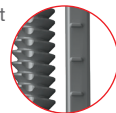

### Ergonomically Advanced Design

Nuance's rectangular design and rubber thumb grips make the syringe feel comfortable to hold and extrude with half the force of leading materials.

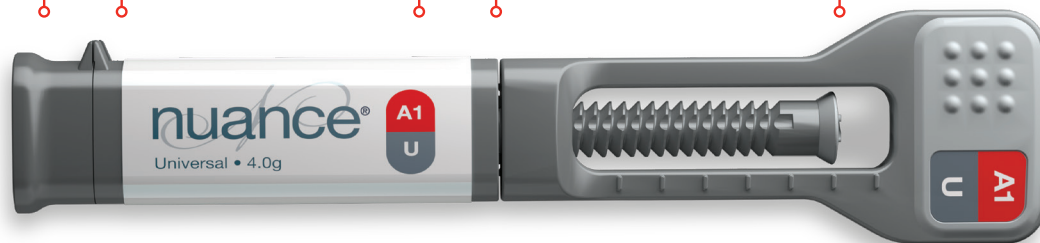

Product	Item #
<b>Nuance Universal Syringe</b>	
Nuance Universal Syringe A1	CR1101
Nuance Universal Syringe A2	CR1102
Nuance Universal Syringe A3	CR1103
Nuance Universal Syringe A3.5	CR1104
Nuance Universal Syringe B1	CR1108
Nuance Universal Syringe B2	CR1109
Nuance Universal Syringe C3	CR1112
Nuance Universal Syringe BL (Bleach White)	CR1106
Nuance Universal Syringe OA2 (Opaque A2)	CR1121
Contents: 1 – 4 g syringe, instructions	

Product	Item #
<b>Nuance Universal Singule</b>	
Nuance Universal Singule A1	CR1135
Nuance Universal Singule A2	CR1136
Nuance Universal Singule A3	CR1137
Nuance Universal Singule A3.5	CR1138
Nuance Universal Singule B1	CR1142
Nuance Universal Singule B2	CR1143
Nuance Universal Singule C3	CR1146
Nuance Universal Singule BL (Bleach White)	CR1140
Nuance Universal Singule OA2 (Opaque A2)	CR1155
Contents: 16 – 0.25 g singules, instructions	

Product	Item #
<b>Nuance Flowable</b>	
Nuance FLOW A1	CR1217
Nuance FLOW A2	CR1215
Nuance FLOW A3	CR1214
Nuance FLOW A3.5	CR1213
Nuance FLOW B1	CR1212
Nuance FLOW B2	CR1211
Nuance FLOW BL1	CR1210
Nuance FLOW BL2	CR1209
Contents: 1 – 1 g syringe, instructions	

Available through



To Order: 1-800-372-4346 8:00am–8:30pm (et)  
[www.henryscheindental.com](http://www.henryscheindental.com)